


BIOLOGICAL MICROSCOPES

N800B – N800T



N800T

Laboratory Biological Microscope

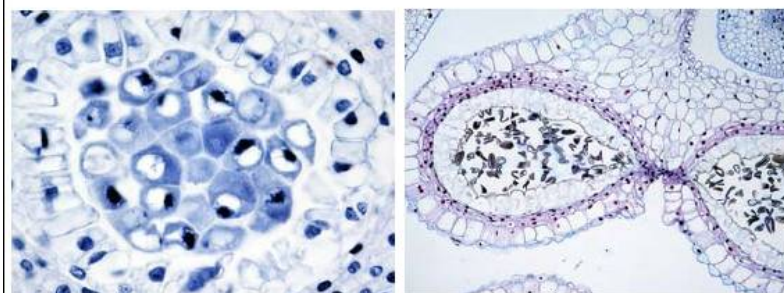

LABOR

Series

SPECIFICATIONS	
Viewing Head	Comfortably Observing – Head inclined at 30° - Suitable for every Observer with Wide Interpupillary Adjustment Range 48-75mm – Diopter Adjustable -5..0..+5.
Eyepiece	Larged Viewing Area with Extra Wide Field Eyepiece EWF 10x/22mm Eyepiece tube Diameter 30mm
Objective	INFINITE PLAN ACHROMATIC OBJECTIVE PL 4x/0.10 PL 10x/0.25 PL 40x/0.65 PL 100x/1.25 OIL
Nosepiece	Big Operatine Space with Backward Quintuple NOSEPIECE with Clickstop.
Condenser	Swing Condenser N.A.: 0.9/0.25 Overcoming Illumination Asymmetry at Low Magnification observation.
Stage	Rectangle Double Layers Mechanical Stage 185x142mm Moving Range 75x55mm – Comfortable Operation with Low position Adjustable Knobs.
Focusing	Coaxial coarse/fine focus system, with tension adjustable and limit stopper the minimum readings is 0.002mm.
Body	Stable and Integral Structure, Comfortable and effective Operation and small Laboratory Space occupancy.
Illuminator	EXTERNAL illumination with Halogen Lamp 12V/30W, Brightness Adjustable
Field Diaphragm	With Aspherical Condenser and Kohler illumination, providing Enough, Even and High Contrast illumination for various Objective Observation.
Alimentazione	110-230 V / 50-60Hz
Accessori in dotazione	Blu Filter 44mm Halogen Lamp 12V 30W Dust Cover

Code	Description
N-800B	Binocular Microscopes PLAN - 30W
N-800T	Trinocular Microscopes PLAN - 30W

Code	Accessories / Description
N811	Eyepiece Pair WF16x
N812	Eyepiece Pair WF20x
N813	Micrometric Eyepiece EWF 10x
N814	Glass Slide 2mm x 0,01mm
N818	Dry Dark field Condenser
N819	OIL Dark field Condenser / Objective 100x for Darkfiled
N820	CCD Adapter 0.5x
N821	Reflex Adapter
N822	Epi-Fluorescent KIT (Torretta – Alimentatore - Lamina)
N823	Objective UIS PL 4x Fluorescent Free
N824	Objective UIS PL 10x Fluorescent Free
N825	Objective UIS PL 20x Fluorescent Free
N826	Objective UIS PL 40x Fluorescent Free
N827	Objective UIS PL100x Fluorescent Free
N828	ULTRAVIOLET Filter System
N829	VIOLET Filter System
N830	Halogen Lamp 6V 30W G4
N831	Temperature Controller Stage



N-800B



ESSE 3 Via Garibaldi 30
14022Castelnuovo D.B. (AT)
tel +39 011 99 27 706
fax +39 011 99 27 506
e-mail esse3@chierinet.it
web : www.4ci.com



Pulse Healthcare Technology
House: 11/1, Shahid Minar Road
Kallyanpur, Dhaka-1216, Bangladesh.
Mobile: +8801708008061
e-mail: info@pht.com.bd
web: www.pht.com.bd

